Acupressuurmat en kussen

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Foreword

About this document

This manual contains all the information for correct, efficient and safe use of the product. Make sure you have fully read and understood the instructions in this manual before installing or using the product. Use the product only as described and only for the purposes specified. Keep this manual for future use.

About us

Feel at home and enjoy optimal comfort with **LifeGoods**. Discover a wide range of quality and affordable products for in and around the home.

From everyday household products to stylish and functional home goods, you'll find everything you need at LifeGoods. Whether it's sturdy party tents, elegant garden chairs, sporty SUP boards or ergonomic office chairs for your (home) workplace, we offer solutions that make your house a home.

Quality, functionality, ease of use and stylish design are our focus. This will make you feel completely comfortable and at home wherever you are. With LifeGoods, you create a home full of comfort and fun, inside and out.

Check out our products at www.lifegoods.nl.

Table of contents

1. Safety	14
2. Introduction	16
2.1 Intended use	16
2.2 Product overview	17
2.3 Product specifications	17
2.4 Package contents	18
3. Instructions for use	18
4. Maintenance and cleaning	18
5. Storage	18
6. Disposal and recycling	18
6.1 Discard	18
7. Service and warranty	19
7.1 Terms of the warranty	19
7.2 Warranty period	19
7.3 Not covered under warranty	20
7.4 Making a claim under the warranty	20
8. Disclaimer	20

1. Safety

WARNING!

Make sure that the packaging material does not get into the hands of children. They can choke on it!

Before use, read the manual and safety information carefully. Follow the instructions and save the manual for the future. The seller is not responsible for failure to follow safety information on the use and handling of the product or damage to property or personal injury.

DANGER

This category contains warnings for situations that could cause serious injury, damage or death if not avoided.

- Do not disassemble, repair or restore the product. This may cause damage to the product. For your own safety, these parts should only be replaced by authorized service personnel. Always send the product to service center for repairs to avoid damage and personal injury.
- Do not use accessories not offered by the supplier. These may pose a danger to the user and damage the product.
- Do not allow children to play with the product. Use and store the product out of the reach of children and pets.

**WARNING**

This category contains warnings for situations that could cause injury or damage if not avoided.

- Note the risks and consequences associated with improper use of the product. Improper use of the product may result in personal injury and damage to the product.
- Use the product only for the purposes indicated in this manual.
- Before each use, check parts and accessories for damage or wear (cracks/breakage) that may adversely affect performance and, if so, replace them with original or manufacturer-recommended parts and accessories.
- If the product is malfunctioning, dropped or damaged, do not use the product and consult our customer service.
- The product is designed and intended for private use. Do not use the product for commercial purposes.
- If you experience pain or discomfort while using the product, stop immediately and consult your doctor.
- Consult a doctor before using the product if you have unexplained pain, are undergoing medical treatment or are using medical equipment.
- Always seek the advice of a professional trainer or doctor to adopt the proper posture for each exercise.

CAUTION

This category contains warnings that need extra attention to avoid situations that could cause minor injuries or damage.

- Keep the product away from sharp and pointed objects that can damage the material.
- Keep the product away from heat, areas where the temperature is relatively high and heat sources such as ovens.
- This product does not serve as a substitute for treatments from a physician or physical therapist. Chronic symptoms may worsen.

2. Introduction


Thank you for choosing LifeGoods product! Always follow safety precautions for safe use.

2.1 Intended use

An acupressure mat is based on the traditional Chinese medicine practice of acupuncture. The mat features plastic points instead of needles, with each point applying pressure to specific points on the body. An acupressure mat can help with relaxation, promote sleep quality and relieve various types of pain, including muscle pain, back pain and headaches. This mat is for domestic use only. It is recommended that you seek medical advice beforehand, especially if you have a medical condition or are pregnant.

The seller is released from liability for any injury or damage caused by the use of the product beyond its intended use. Any modification of the product may affect your safety and warranty.

2.2 Product overview

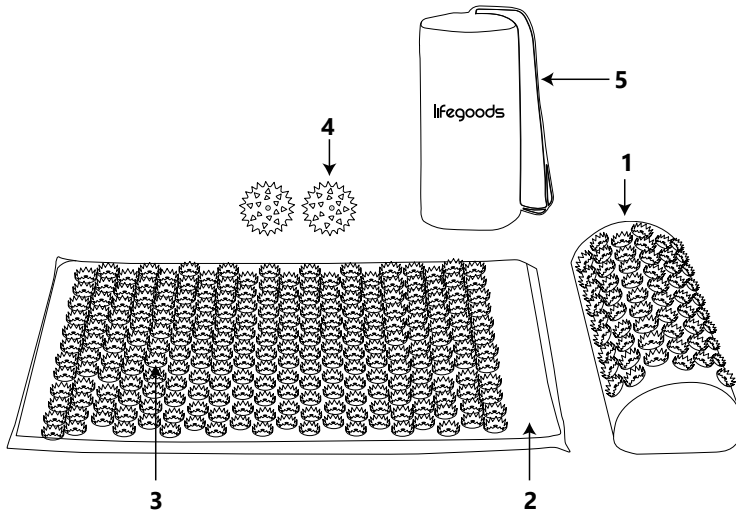1. Pillow

2. Mat

3. Acupressure points

4. Massage ball

5. Carrying bag

2.3 Product specifications

Dimension mat	66 x 40 x 2 cm
Pillow size	37 x 15 x 10 cm
Dimension of carrying bag	58 x 31.5 cm
Weight	0.6 kg
Material	Polystrene, polyester and cotton
Number of points mat	230 x 27
Number of points pillow	66 x 27

2.4 Package contents

Check that the package contents match the list in the manual. Notify the seller if parts are missing. If the product appears damaged, immediately file a complaint with the carrier and provide the seller with a detailed description of the damage. Keep the damaged parts and packaging to substantiate your report.

- 1x Acupressure mat and pillow
- 1x Carrying bag
- 1x Instruction booklet
- 1x Pillow
- 2x Massage ball

3. Instructions for use

See attached instruction booklet.

4. Maintenance and cleaning

The product is not suitable for the washing machine. Clean the product with a damp, clean cloth and let it air dry. Do not use abrasive or aggressive cleaning agents.

5. Storage

Store the product in a dry, clean place. Do not store the product in extremely high or low temperatures.

6. Disposal and recycling

Think about the environment and contribute to a cleaner living environment! You can contact your municipality for all information on disposal options for disposed products.

6.1 Discard

At the end of its life, do not dispose of the product with normal household waste, but hand it in at a government designated collection point for recycling.

Dispose of the packaging and accessories in accordance with the environmental regulations in force in your municipality. Observe the symbols on the various packaging materials and dispose of them separately – if necessary. Packaging materials are labeled with letters and numbers: 1-7: plastic, 20-22: paper and cardboard, 80-98 composite materials.

Separate waste into glass, paper and plastic for recycling with Tri sélectif. Reduce landfill waste and conserve natural resources. Contribute to sustainability by consciously separating.

7. Service and warranty

If you need service or have a question about your product, please contact our customer service team at service@lifegoods.nl.

You get a 2-year warranty on this product from the date you purchased it. If the product breaks down, you are always entitled to legal protection through the seller. This warranty does not change those legal rights.

7.1 Terms of the warranty

The warranty starts on the day you bought the product. Therefore, keep the original receipt or email confirmation carefully – this is your proof of purchase. If the product breaks down within two years of purchase due to a material or manufacturing defect, we will repair or replace the product free of charge. For this you must:

- Contact our customer service department at service@lifegoods.nl.
- Enclose a copy of the proof of purchase.
- Briefly write down what the problem is and when it arose.

If the defect falls under warranty, you will receive the repaired or a new product back. This does not extend or restart the warranty period.

7.2 Warranty period

The warranty lasts for 2 years. This period is not extended if your product is repaired or replaced. Are there any damages or defects on the product when it is unpacked? Please report this to our customer service department immediately. Repair or replacement outside the warranty period are not free of charge.

7.3 Not covered under warranty

The product is carefully made and tested. The warranty applies only to errors in material or manufacturing. The following are not covered by the warranty:

- Parts that wear out from normal use.
- Fragile parts getting damaged.
- Damage caused by misuse, improper maintenance or if you have used the product in a manner not specified in the instructions for use.
- Use for business or commercial purposes.
- Damage caused by abuse, violence or modification by anyone other than our own service department.

Read the manual carefully and follow the instructions.

7.4 Making a claim under the warranty

Follow these steps for quick processing:

- Have the proof of purchase, item number and EAN code at hand.
- You can find the part number on the product nameplate, in the user manual (at the bottom of the title page), or on a sticker on the back or bottom of the device.
- If you have any problems, please contact us first at service@lifegoods.nl or via the contact form at www.lifegoods.nl.
- After consulting with our customer service department, you can send the product to us if necessary. Please use the return address provided by us on the return label. Include a copy of the proof of purchase and write what the defect is and when it occurred.

Important

- Ship the product only to the address provided.
- Make sure the packaging is sturdy.
- Do not use express shipping, freight shipping or special forms of delivery.
- Also send all accessories that came with your purchase.

8. Disclaimer

Changes subject to change; specifications are subject to change without notice.

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